



mediroyal 

Orthoses are our passion

SRX  [®]

BY MEDIROYAL

NEW!

*Support during Activity,
taking you further.*

SRX[®] is a registered trademark
of Mediroyal Nordic AB, Sweden

© Mediroyal Nordic AB, Sweden 2021






BY MEDIROYAL

The SRX® line has been designed to support people in activity. The new innovative SRX® material has a 3-dimensional top surface textile that provides both elasticity and compression. The core is made from CR-foam that contains closed micro air-cells, to provide a lightweight material with good heat insulation and no moisture absorption.

Our Neptune™ material, which is hydro-dynamic, comfortable and quick-drying, has been laminated on the inside. The new materials, in combination with a functional design, and our Scandinavian heritage that values an active lifestyle, are all essential parts of our SRX® line.

Support during Activity, taking you further.

mediroyal

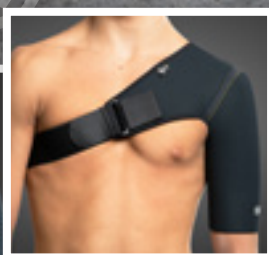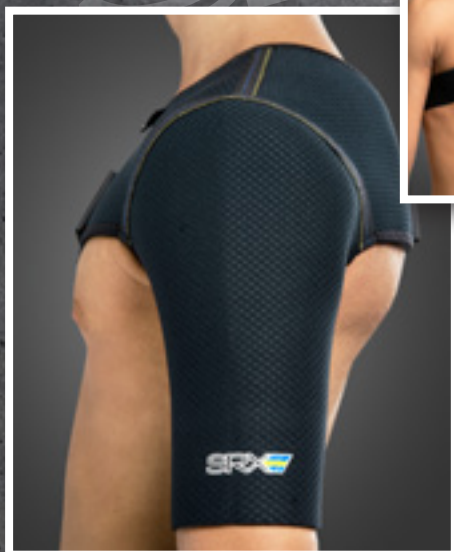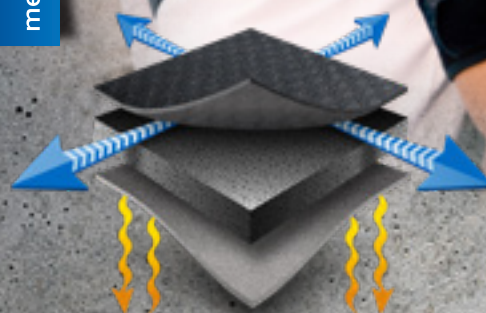

SRX810 · SRX® Shoulder Support

This shoulder support has been designed with the 3 mm SRX® material that provides support, compression, and local heat insulation. The design makes it low profile, good fit and very functional. The closure can be alternated left or right by easily changing the attachment. The shoulder support can be used for local inflammations, mild instability, or temporary pain caused by muscle strains.

Size	XS	S	M	L	XL
cm	24–28	28–32	32–35	35–38	38–41

The circumference measurement is taken around the center of the upper arm.



SRX811 · SRX® Shoulder Support X

This shoulder support is made from a 3 mm SRX® material that is combined with an individual cross-strap over the shoulder joint to provide a more individual support and stability. The upper arm has an NRX® material with a friction effect along with a hook adjustment. The support can be worn on both left and right shoulder by alternating the attachment. The shoulder support can be used for local inflammations, mild instability, or temporary pain caused by muscle strains.

Size	XS	S	M	L	XL
cm	24–28	28–32	32–35	35–38	38–41

The circumference measurement is taken around the center of the upper arm.



SRX820 · SRX® Elbow Support

This elbow support has been designed with the 3 mm SRX® material to provide support, compression, and local heat insulation over the elbow joint. The functional design, made from 3 separate pieces and a plain piece over the elbow crease, provides a low profile. The support has hook and loop closure over the upper and lower part for more adjustability. The elbow support can be used for local inflammations, muscle strains and mild instabilities.

Size	XS	S	M	L	XL
cm	20–23	23–25	25–28	28–31	31–33

The circumference measurement is taken around the center of elbow joint.



SRX821 · SRX® Universal Elbow Support

This elbow support is made from a 3 mm SRX® material and has a universal wrap-around design that is applied over the elbow joint. It provides support, compression, and local heat insulation. The hook and loop closure makes it possible to off-set a larger circumference in the upper arm with a smaller circumference over the lower arm. The top surface is made from a loop fabric that allows attachment of the hook parts anywhere. The elbow support can be used for local inflammations, muscle strains and mild instabilities.

Size	Universal
cm	20–30

The circumference measurement is taken around the center of elbow joint.



SRX822 · SRX® Universal Elbow Sleeve

This elastic elbow sleeve is made from a circular knitted fabric with spiral elastic yarn. It provides a light support and compression over the elbow joint. The elbow sleeve can be used for local inflammations, muscle strains and mild instabilities.

Size	Universal
cm	22–28

The circumference measurement is taken around the center of elbow joint.





SRX823 · SRX® Universal Tennis Elbow Strap

The tennis elbow strap is made from a 3 mm SRX® material with Neptune™ lining on the inside. It has a fixed pad just under the D-ring attachment to provide a light support over the muscles. The strap is elastic which provides a dynamic compression to hold it in place during activity. The hook and loop closure makes it versatile and flexible in adjustment. The strap can be used for both tennis elbow as well as golfers elbow problems.

Size	Universal
cm	20–32

The circumference measurement is taken 5 cm below the elbow joint.



SRX215 · SRX® Universal Firm Tennis Elbow Strap

The firm tennis elbow strap is made from a double laminated loop fabric that is non-elastic. It has a pad that is attached with hook on the inside that allows individual positioning. The firm strap provides a more distinct support over the muscles and the hook and loop closure makes it versatile in adjustment. The strap can be used for both tennis elbow as well as golfers elbow problems.

Size	Universal
cm	20–32

The circumference measurement is taken 5 cm below the elbow joint.



SRX824 · SRX® Universal Wrist Strap

This wrist strap is made from 3 mm SRX® material with Neptune™ lining on the inside. The strap around the wrist is elastic for a good support and compression and the hook and loop closure makes it versatile and flexible in adjustment. The wrist strap can be used to support the wrist for local strains and instabilities as well as for tennis elbow problems, to reduce wrist movement and off-load the muscles in the lower arm.

Size	Universal
cm	12–20

The circumference measurement is taken around the wrist.



SRX825 · SRX® Universal Wrist & Thumb Strap

This wrist & thumb strap is made from 1,5 mm SRX® material with loop surface on the top and Neptune™ lining on the inside. The thumb loop is attached over the thumb and the strap is applied around the wrist, providing a light support for both the thumb and wrist. The hook and loop closure makes it versatile and flexible in adjustment. The wrist strap can be used for local inflammation, overloading, strains and mild instabilities in the wrist and thumb.

Size	Universal
cm	12–20

The circumference measurement is taken around the wrist.



SRX826 · SRX® Universal Wrist Strap

This wrist wrap is made from 3 mm SRX® material with a loop surface on the top and Neptune™ lining on the inside. The design features a thumb hole for easy application over the thumb, then the wrap is applied around the wrist and the elastic strap attached as a reinforcement. The hook and loop closure makes it both stable, versatile, and yet flexible in adjustment. The wrist strap can be used to support the wrist for local strains and instabilities as well as for tennis elbow problems, to reduce wrist movement and off-load the muscles in the lower arm.

Size	Universal
cm	12–22

The circumference measurement is taken around the wrist.





BY MEDIROYAL


Support during Activity, taking you further.

mediroyal




SRX840 · SRX® Back Support

This back support is made from 5 mm SRX® material that provides support, and local heat insulation to the lower back area. The inside is lined with the Neptune™ material for high comfort and quick drying. The hook and loop closure in the front makes it easy to apply and adjust. The back support can be used for temporary low back pain, muscle strains, or to provide a light support.

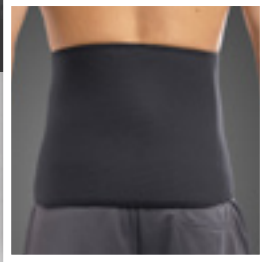

Size	S	M	L	XL	XXL
cm	70–80	80–90	90–100	100–110	110–120

The circumference measurement is taken around the waist.

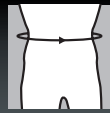



SRX841 · SRX® Back Support Stable

This back support is made from a ventilated elastic with a low profile to be able to use it under clothes. The back part has plastic stays to provide support for the lower back. The inside of the lower back part has an anti-slip function to prevent migration of the back support. As an additional reinforcement there is a double strap system attached in the back, that can be fastened in the front to provide individual support and compression. The back support can be used for temporary low back pain, muscle strains, or to provide support to the lower back during activity.

Size	S	M	L	XL
cm	60–75	75–90	90–120	120–150

The circumference measurement is taken around the waist.

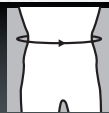

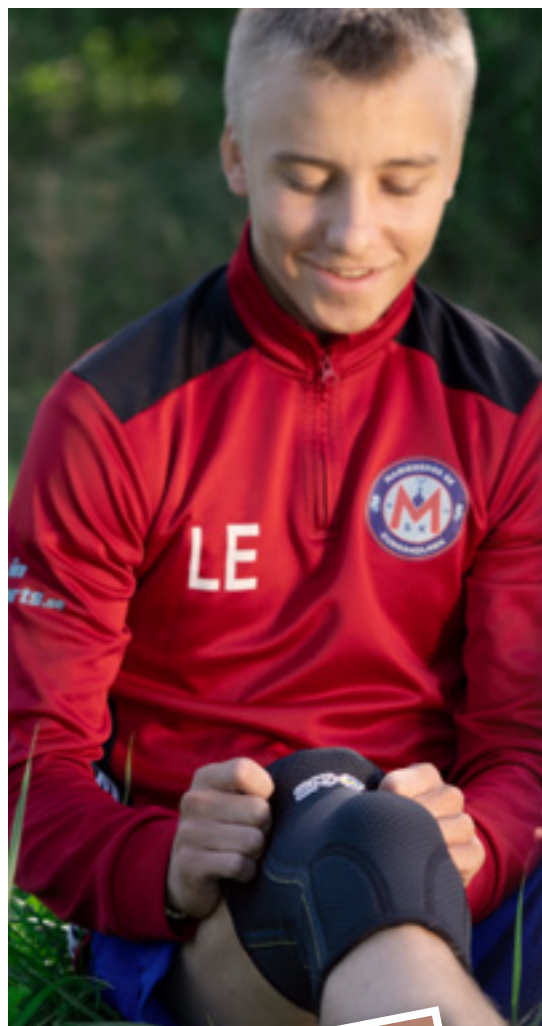

SRX852 · SRX® Multi Sport Knee

This model is made from a combination of 5 mm SRX® material in the main body and a 3 mm SRX® material in the knee fold. The anatomical design provides a good support, very good fit and the side panels makes the knee support more flexible during activity, especially for multi-sport or cross fit. The SRX® material is lined with Neptune™ which provides a high comfort and dries quickly. The knee support can be used for local inflammations, overuse, local swelling, mild instabilities, or temporary knee pain.

Size	S	M	L	XL
cm	33–36	36–39	39–42	42–45

The circumference measurement is taken around the knee joint.



SRX860 · SRX® Impact Knee Support

This model has been designed to protect the knee joint from impact during sports by using a pad made from shock absorbing CR15-foam, laminated with Kevlar®. The basic model is the same as SRX852 but with a slightly shorter design in the front. The anatomical design provides a good support, very good fit and the side panels makes the knee support more flexible during activity. The Neptune™ material is laminated on the inside which provides high comfort and dries quickly. The knee support can be used for local inflammations, overuse, local swelling, mild instabilities, temporary knee pain or to protect the knee joint from impact.

Size	S	M	L	XL
cm	33–36	36–39	39–42	42–45

The circumference measurement is taken around the knee joint.



SRX854 · SRX® Knee Strap

This model is made from a 3 mm SRX® material with an anti-friction surface on the inside and a loop fabric on the top. The anti-friction surface provides a good skin grip and prevents the strap from migrating during activity. The strap has a pad in the front that is placed over the kneecap tendon to provide support. There is an additional hook strap that secures the strap and provides individual support to the pad. The knee strap can be used for local inflammations, jumpers knee problems as well as Osgood–Schlatter disease.

Size	XXS	XS	S	M	L
cm	26–29	29–32	32–35	35–38	38–42

The circumference measurement is taken around the knee joint.



SRX855 · SRX® Universal Knee Strap

This model is a simple, versatile, and universal knee strap that is made from SRX® material with a top surface made from loop fabric. The design makes it ideal as an entry strap for kneecap tendon problems. The hook and loop adjustment makes it versatile and easy to apply. The pad is placed over the kneecap tendon providing a light support. The compression is adjusted by the circumference tension of the strap. The knee strap can be used for local inflammations, mild jumpers knee problems as well as mild Osgood–Schlatter disease.

Size	Universal
cm	29–42

The circumference measurement is taken around the knee joint.




SRX850 · SRX® Knee Support

This knee support is made from a 3 mm SRX® material lined with Neptune™ on the inside for optimal comfort. The anatomical design, that is slightly flexed, provides a good fit and function as well as local heat insulation during activity. The knee support can be used for local inflammations, overuse, local swelling, mild instabilities, or temporary knee pain.

Size	XS	S	M	L	XL	XXL
cm	29–32	32–35	35–38	38–42	42–46	46–50

The circumference measurement is taken around the knee joint.



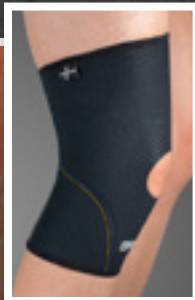

SRX851 · SRX® Knee Support Open Kneecap

This knee support is made from a 3 mm SRX® material lined with Neptune™ on the inside. The design has an opening for the kneecap to provide pressure off-loading. The anatomical design provides a good fit and function as well as local heat insulation during activity. The knee support can be used for local inflammations, overuse, local swelling, mild instabilities, or temporary knee pain.

Size	XS	S	M	L	XL	XXL
cm	29–32	32–35	35–38	38–42	42–46	46–50

The circumference measurement is taken around the knee joint.



SRX856 SRX® Knee Support Wrap

This model is a shorter wrap-round design that is very versatile and can accommodate an over-sized thigh or calf. The model is made from 3 mm SRX® material with a loop fabric surface and Neptune™ lining on the inside. To provide support to the knee there are two spiral stays on the side of the knee along with two elastic, detachable straps that can be positioned anywhere on the surface. The knee support can be used for local inflammations, overuse, local swelling, mild instabilities, or temporary knee pain.

Size	1	2
cm	29–38	38–52

The circumference measurement is taken around the knee joint.



SRX857 · SRX® Knee Support Spiral

This circular shaped model is made from 3 mm SRX® material lined with Neptune™ on the inside. The model has two spiral stays on each side of the knee joint to provide additional support and an opening for the kneecap. Two elastic and detachable straps provide more support and individual compression over the knee joint. The knee support can be used for local inflammations, overuse, local swelling, mild instabilities, or temporary knee pain.

Size	XS	S	M	L	XL	XXL
cm	29–32	32–35	35–38	38–42	42–46	46–50

The circumference measurement is taken around the knee joint.



SRX858 · SRX® Hinged Knee Brace

This circular shaped model with an opening for the kneecap is made from 3 mm SRX® material lined with Neptune™ on the inside. The model has two lightweight hinges made from composite and aluminium to provide support and stability for the knee. Two elastic and detachable straps can provide more individual compression over the knee joint. The knee brace can be used for local inflammations, overuse, local swelling, or moderate instabilities caused by ligament or meniscus injuries.



Size	XS	S	M	L	XL	XXL
cm	29–32	32–35	35–38	38–42	42–46	46–50

The circumference measurement is taken around the knee joint.



SRX859 · SRX® Hinged Knee Wrap Brace

This model is a wrap-around design that is very versatile and can accommodate an over-sized thigh or calf. The model is made from 3 mm SRX® material with a loop fabric surface and Neptune™ lining on the inside. The model has two lightweight hinges made from composite and aluminium to provide support and stability for the knee. Two elastic and detachable straps can provide more individual compression over the knee joint. The knee brace can be used for local inflammations, overuse, local swelling, or moderate instabilities caused by ligament or meniscus injuries.

Size	1	2
cm	29–38	38–52

The circumference measurement is taken around the knee joint.





SRX870 · SRX® Calf Support

This calf support is made from 1,5 mm SRX® material with Neptune™ lining on the inside. The material has a high stretch with controlled compression providing good support to the calf muscle. The curved anatomical shape over the calf makes it possible to accommodate muscles that increase in volume during activity. The compression in the support is different depending on the size selected. A size smaller will provide higher compression. The lower part of the support has three stop-seams that makes it possible to trim the length of the support if needed. The calf support can be used for pain or muscle strains, local inflammations, shin-splints, and Achilles' problems.


Size	XS	S	M	L	XL	XXL
cm	33–34	34–37	37–40	40–42	41–44	43–46

The circumference measurement is taken around the widest part of the calf. Note that the shape of the calf muscle can influence the individual fit of the support. The length of the support can be trimmed over the Achilles if it is too long. Always cut below the stop-seam.



SRX871 · SRX® Ankle Support

This ankle support is made from 3 mm SRX® material with Neptune™ lining on the inside. The material provides good support and compression in combination with local heat insulation. An elastic strap is attached under the foot which can be used to provide additional support and compression over the ankle joint. The ankle support can be used for local heat insulation, light support, and compression over the ankle joint after mild sprains.


Size	XS	S	M	L	XL	
cm	19–20	20–21	21–23	23–25	25–27	
Shoe size EUR	35–37	38–40	40–42	42–44	44–46	

The circumference measurement is taken around the ankle joint, over the Malleolus. Shoe size can also be used as reference. If the measurements do not match up, we recommend using the circumference measurement.



SRX872 · SRX® Universal Ankle Support


This universal ankle support is made from a circular knitted fabric with spiral elastic yarn. It provides a light support and compression. An elastic strap that is attached under the foot which can be used to provide additional support and compression over the joint. The ankle support can be used for light support and compression over the ankle joint after mild sprains.

Size	Universal	
cm	18–24	
The circumference measurement is taken around the ankle joint, over the Malleolus.		



SRX875 · SRX® Ankle Brace Stable

This stable lace-up ankle brace is made from a ballistic nylon material, lined with an air-mesh material. Over the Achilles there is an elastic stretch zone to prevent pressure on the tendon. The tongue is made from SRX® pressure absorbing CR-foam. The nylon straps attached on the side of the brace should be crisscrossed under the foot and attached to each side of the brace to provide individual stability. The upper wide strap secures both the straps and lace in place once adjusted. The ankle brace can be used for mechanical support and stability over the ankle joint after mild to moderate sprains or to provide support for an instable ankle during activity.

Size	XS	S	M	L	XL	
cm	19–20	20–21	21–23	23–25	25–27	
Shoe size EUR	35–37	38–40	40–42	42–44	44–46	

The circumference measurement is taken around the ankle joint, over the Malleolus. Shoe size can also be used as reference. If the measurements do not match up, we recommend using the circumference measurement.





SRX880 · SRX® Thermal & Compression Pants

The SRX® Thermal & Compression pants is made from a combination of a 0,7 mm ultra-stretch CR-material and 1 mm SRX® compression material, both lined with Neptune™ on the inside. The 1 mm SRX® compression material is made with a higher density and is placed over the front side of the thigh, running up over the groin area in the front to provide higher compression and support. The back part of the pants is slightly higher to provide more heat to the lower back area. An adjustable drawstring is positioned over the waist. The pants can be used to provide local heat insulation and compression over the buttocks, thighs and groin which can improve circulation and elasticity of muscles.

Size	XS	S	M	L	XL
cm	60–75	70–85	80–95	90–105	105–115

The waist size measurement is an indication for the correct size however the body constitution can influence the fit. The material has an ultra-stretch which means that the same size could overlap the sizing.



SRX990 · SRX® Ice-Wrap for Knee, Ankle, Elbow and Shoulder

The SRX® Ice-Wrap is a multi-function wrap with an icebag that can be filled with ice and water to provide cold to an injured body part. The product contains an icebag and a wrap with two long straps that will hold it in position, making it possible to apply over shoulder, elbow, knee, or ankle. The Ice-Wrap should only be used for a limited time of 20-30 minutes and never directly on bare skin. Always use a towel to protect the skin. The ice-wrap can be used after an acute injury to administrate cold therapy for up to 30 minutes to the shoulder, elbow, knee, and ankle joint.

Universal size



SRX700 · SRX® Compression Socks · White SRX701 · SRX® Compression Socks · Black



The SRX® Sport Compression Socks have been developed together with the Swedish Handball Federation for the Swedish Men's and Women's National Handball teams.

The sock has a graded medical compression, equivalent to compression class 2, to battle muscle fatigue in the calf muscles during endurance sports. They are made from a unique fiber composition that combines durable and comfortable bamboo fibers with an active antimicrobial effect from charcoal, together with Nylon yarn. The elastic yarn used is latex-free Lycra. The SRX® sock has an anatomical design with a larger calf area and a reinforced ankle part, where the compression grades off from. Both the toe and the heel parts have been reinforced for a good durability during extreme sports. The upper part of the sock has a flexible zone to prevent discomfort where the sock ends. The length of the sock has been made flexible to fit both normal shins as well as long ones.

The compression socks can provide less fatigue muscles, improve healing, prevent edema after trauma and enhance the range of motion in the joint.



Size	Shoe size	Calf measurement	Ankle measurement
XS	34–37	23–33 cm	16,5–19 cm
S	37–39	28–38 cm	19–21,5 cm
M	40–43	33–45,5 cm	21,5–24 cm
L	44–47	35,5–46 cm	23–26,5 cm
XL	48–	38,5–51,5 cm	26,5–32 cm



Quick drying



Speeds up recovery



Good insulation ability



Increases circulation



Antibacterial & natural deodorizer



UV resistant

SRX042 · SRX® Sport Bag · Black

A water repellent 42 liter sport bag that can be used as both a regular bag and a back-pack. Made from high quality ballistic nylon and durable TPU. The shoulder straps have quick-buckle release that can be removed easily and stored in the top lid on the bag. Transparent pocket on the top for a business card. SRX® by Mediroyal print on one of the sides of the bag.





Mediroyal Nordic AB is a Swedish company focused on medical products, with more than 20 years' experience in developing functional orthotic solutions.

We produce our products and materials in our own production facilities in Taiwan, Republic of China and export to over 50 countries around the world.

Orthoses are our passion

mediroyal 

Mediroyal Nordic AB
Box 7052. Staffans väg 6B · SE-192 07 Sollentuna. Sweden
info@mediroyal.se · www.mediroyal.se

Mediroyal Asia Pacific Co. Ltd.
Unit A10, 20F., No. 282, Shizheng N. 2nd Rd., Xitun Dist., Taichung City 407, Taiwan (R.O.C.)
+886-4-2254-8069 · sales.mediroyal@gmail.com · www.mediroyal.se/en