




mediroyal 

Orthoses are our passion

Royal Putty Nova & Featherlite Nova

We have further developed our popular hand therapy putty after suggestions from users and our customers. The two new Nova-putties have a drier feeling and are less oily compared to our previous qualities and the dynamic resistance has also been enhanced.


**FURTHER
DEVELOPED
FORMULA**

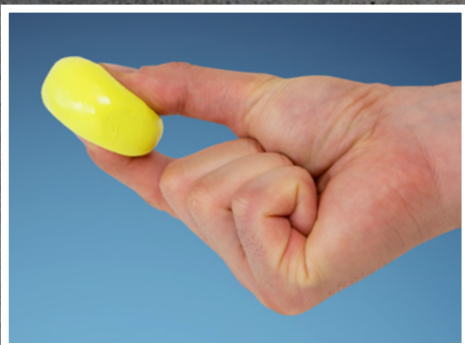

The composition of the putty is silicone elastomers with different densities. In our new qualities we have increased the dynamic resistance by cross-linking the elastomers. When training with a higher intensity the putty responds with a higher resistance. At low intensity the putty responds with a lower resistance. This makes it possible to use the same resistance level longer during the rehabilitation phase.

The softer resistance levels have been modified to work better for neurological and rheumatic patients.

We are also introducing a middle package of 450 g/335 cc as a complement to our existing packages.

See resistance variants
and package sizes 

© Mediroyal Nordic AB 2022

Mediroyal Nordic AB · Box 7052 · 192 07 Sollentuna · Sweden · Tel +46 8 506 766 00
info@mediroyal.se · facebook.com/Mediroyalnordic · www.mediroyal.se



Royal Putty Nova

Available in six resistance levels. The lower resistance levels, XX-Soft and X-Soft have been modified to work better for neurological and rheumatic patients. Even with a low resistance it does not smudge.

All the resistance levels can be blended for an individual resistance.

Resistance	Color	85 g	450 g	2,25 kg
XX-Soft	White	1111	1131	1121
X-Soft	Beige	1112	1132	1122
Soft	Yellow	1113	1133	1123
Medium	Red	1114	1134	1124
Firm	Green	1115	1135	1125
X-Firm	Blue	1116	1136	1126



Featherlite Nova

The new enhanced Featherlite Nova is much lighter than the previous quality. The dynamic resistance has been improved as well as the elasticity to make it more functional for hand rehabilitation.

Featherlite Nova is not sold in grams but in volume cc due to the low weight.

Resistance	Color	65 cc	335 cc	1625 cc
X-Soft	White	2111	2131	2121
Soft	Yellow	2112	2132	2122
Soft-Plus	Lime	2113	2133	2123
Medium	Green	2114	2134	2124
Firm	Blue	2115	2135	2125

