



AFO Dorsal Splint

AFO Dorsal Splint is a dorsal splint in polypropylene plastic. Supplied with pre-assembled Velcro strap with calf padding and adjustable ankle strap. Outlets in the heel area provide a better fit and make the orthosis work well in the shoe. Can be shaped and customized using heat for a better fit. Places the foot in a neutral position and prevents drop foot.



MR5801 • AFO Dorsal Splint

- PP splint with open heel
- Prevents drop foot
- Lifts and supports the foot
- Flexible design
- Fits well in most types of shoes
- Thin foot part that is easy to cut with scissors, can be individually adjusted

Shoe size • Available in left and right versions

S = 34-37

M = 38-41

L = 42-45



AFO Dorsal Splint


-without band



MR5800 • AFO Dorsal Splint without band

Foot drop splint without application tape can be adapted perfectly with several of our materials that are available on a roll. In the picture you can see how we have used Neoprene with velcro textile for padding at the back. For easy and quick application without having to fasten the band with rivets, we use the NRX Low Profile Hook which is self-adhesive. Around both ankle and calf, soft velcro, low profile is used.

Sizes: see previous page

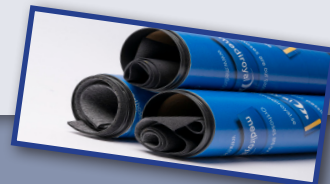

[EU490101](#) • Neoprene with OP loop textile/Black Nylon • 1 mm • W: 58 cm • L: 100 cm
[EU445060](#) • NRX Low Profile Hook, self Adhesife • Transparent • 50 mm x 25 m
[EU445062](#) • Low profile, Soft Velcro • 50 mm x 25 m