

Why choose Go up BASIC

- 1 Perfect fit on the patient
- 2 Light-weight and easy to handle
- 3 New "PUSH&SET" range of motion control for easy adjustment

Why choose Go up PRO

- 1 Telescopic
- 2 Easy to adapt to the patient
- 3 New "PUSH&SET" range of motion control for easy adjustment

9000P0019TEN001.0124

go up BASIC

REF. 24045
Post-operative knee orthosis with graduated movement

Code	REF.24045
Size	universal
Length of the hinges cm	56
Color	grey

ambidextrous



CUSTOMIZABLE FIT



EASY ADJUSTING



SINGLE PATIENT MULTIPLE USE

go up PRO

REF. 24046
Telescopic post-operative knee orthosis with graduated movement

Code	REF.24046
Size	universal
Length of the hinges cm	50-65
Color	grey

ambidextrous



CUSTOMIZABLE FIT



EASY ADJUSTING



SINGLE PATIENT MULTIPLE USE



TELESCOPIC

Indications

- Post-operative treatment requiring the protection of the knee's flexion-extension (ligament reconstruction, anterior tibial apophysis transposition, fractures of the distal third of the femur, fractures of the proximal third of the tibia, patella fractures, tendon injuries of the extensor apparatus)
- Post-traumatic phase of knee sprains.
- Varising or valgus osteotomies



orthoservice.com



Headquarter: ORTHOSERVICE AG
Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland
Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89
info@orthoservice.com · www.orthoservice.com

ORTHOSERVICE DEUTSCHLAND GmbH
Flugstraße 8 · D-76532 Baden-Baden · Deutschland
Tel. 0049 (0) 7221 991 39 11 · Fax 0049 (0) 7221 991 39 13
info@orthoservice.de · www.orthoservice.de

ORTHOSERVICE POLSKA Sp. z o.o.
ul. Warszawska 416a · 42-209 Częstochowa · Polska
Tel. +48 (0) 34 340 13 10 · (NIP): 9492246785
info@orthoservice.pl · www.orthoservice.pl

RO+TEN S.r.l.
Sede legale: Via Marco De Marchi 7 · I-20121 Milano (MI) · Italia
Sede operativa e amministrativa:
Via Comasina 111 · I-20843 Verano Brianza (MB) · Italia
Tel. 0039 039 601 40 94 · Fax 0039 039 601 42 34
info@roplusten.com · www.roplusten.com



Go up BASIC

Go up PRO

go up BASIC / PRO

Post-operative
knee orthosis with
graduated movement



ORTHOSERVICE
RO+TEN

Take care feel better




ORTHOSERVICE
RO+TEN

Take care feel better

New range of motion control "PUSH&SET"

Enables an intuitive and safe flexion-extension adjustment

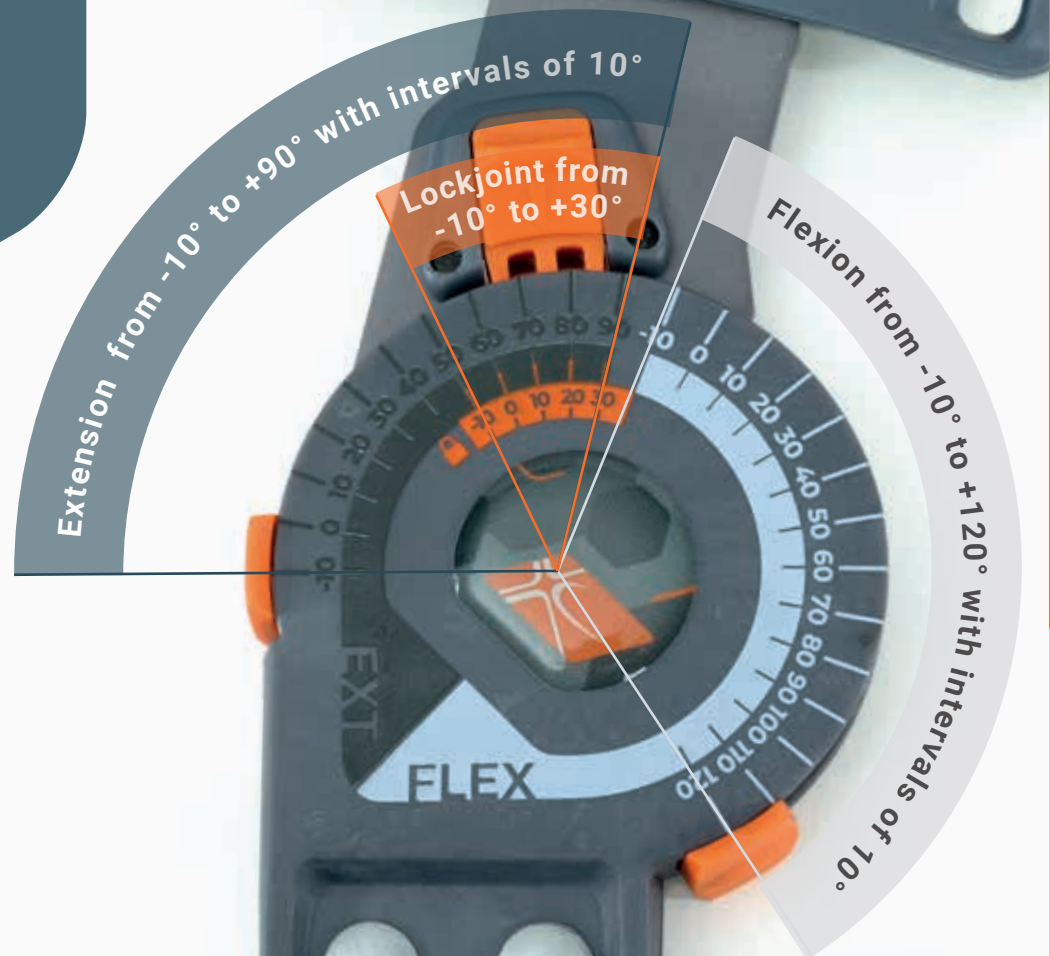

Extension from -10° to +90° with intervals of 10°

Lockjoint from -10° to +30°

Flexion from -10° to +120° with intervals of 10°



Specific button for fast lock joint from -10° to +30° with intervals of 10°

go up BASIC



LATERAL ARTICULATED HINGE IN A LIGHT ALUMINIUM ALLOY
for optimal maneuverability

REMOVABLE AND WASHABLE CONDYLE PROTECTION PADS

SINGLE PIECE MTP PADS
to limit brace movement and protect the leg from direct contact with the hinges

UNIVERSAL SIZE
Enables to work with minimum stocks in your warehouse



UNIVERSAL CIRCUMFERENCE
for a perfect fit on the patient

go up PRO



TELESCOPIC AND EASY TO ADJUST ON THE PATIENT
thanks to the buttons present on the valves

HINGES IN A PLASTIC MATERIAL
with printed Velcro®

STRAP LOOP VALVES
with structural lock joint

REMOVABLE AND WASHABLE CONDYLE PROTECTION PADS

UNIVERSAL SIZE
hinge length can be regulated



MTP PADS IN STRIPS
to limit brace movement and protect the leg from direct contact with the hinges

RENEWED DESIGN

Two new post-operative at your service