

# The all-time best support for the knee

Designed to support knee instability or torn ligaments, it represents a first choice brace for **any instability of the knee**, whether it is due to traumatic or degenerative reasons. Recommended for both **preop and postop phases**, as a conservative treatment or as **preventive protection**.

## Front fit

The Pluspoint Fast knee brace is put on from the front without having to insert the foot. This feature makes it extremely **practical and easy to use**.

## THE STRAPS ARE INCREDIBLY EASY TO FASTEN

The straps don't need to be removed and reinserted every time, because they fasten with **small levers**, which maintain the regulation made on the first application.



## EFFICIENCY OF THE 4 THRUST POINTS

The 4 thrust points system, guarantees efficient support, providing optimal stabilization of the knee.

9DCDP01057E/A001.0322

# pluspoint FAST

## Functional Knee orthosis

REF. 22040

### Indications

- ACL, MCL, LCL, PCL instability and ruptures, combined instability
- Complex instabilities of the knee joint, traumatic or degenerative
- Pre/Post-surgery both rehabilitative or conservative
- Ligament protection during sports



| Item                                      | REF.22040 |       |       |       |       |
|---|-----------|-------|-------|-------|-------|
| Size                                      | S         | M     | L     | XL    | XXL   |
| Thigh circ. taken 15 cm above mid-patella | 39/47     | 47/53 | 53/60 | 60/67 | 67/75 |
| Length cm                                 | 40        | 40    | 40    | 40    | 40    |
| Colour                                    | dark grey |       |       |       |       |

Specify right or left



orthoservice.com



FSC® C003498



UN EN ISO 13485:2016  
UN EN ISO 9001:2015



**Headquarter: ORTHOSERVICE AG**  
Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland  
Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89  
info@orthoservice.com · www.orthoservice.com  
**Niederlassung Deutschland: Orthoservice Deutschland GmbH**  
Flugstraße 8 · D-76532 Baden-Baden · Deutschland  
Tel. 0049 (0) 7221 991 39 11 · Fax 0049 (0) 7221 991 39 13  
info@orthoservice.de · www.orthoservice.de  
**Siedziba w Polsce: Orthoservice Polska Sp. z o. o.**  
ul. Warszawska 416a · 42-209 Częstochowa · Polska  
Tel. +48 (0) 34 340 13 10  
info@orthoservice.pl · www.orthoservice.pl  
**Sede italiana: RO+TEN S.r.l.**  
**Sede legale:** Via Marco De Marchi, 7 · I-20121 Milano (MI) · Italia  
**Sede operativa e amministrativa:**  
Via Comasina, 111 · I-20843 Verano Brianza (MB) · Italia  
Tel. 0039 039 601 40 94 · Fax 0039 039 601 42 34  
info@roplusten.com · www.roplusten.com

Società soggetta a Direzione e Coordinamento (art. 2497bis CC):  
Orthoservice AG (CH) · 6830 Chiasso (TI) · Switzerland

# pluspoint FAST

## Functional Knee orthosis



ORTHOSERVICE  
RO+TEN

Take care feel better




ORTHOSERVICE  
RO+TEN

Take care feel better

# Shapable anatomical frame with front fit

Made of a **special aluminum alloy**, it is shapable, especially on the shine grip and thigh. The **front fit** makes it easy for the patient: fast and intuitive to wear.

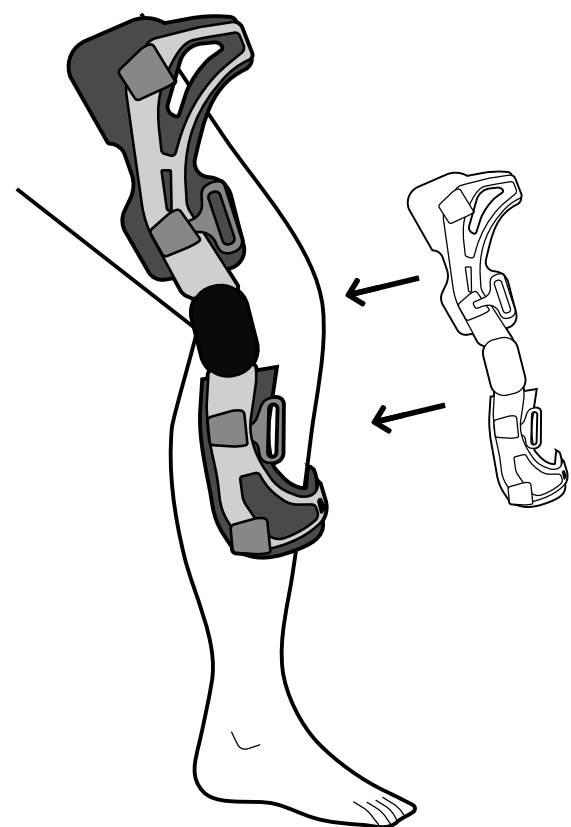

### CLOSURE WITH LEVER BUCKLES

The straps are not inserted in the buckles but are hooked to the frame by means of **special lever buckles**, also equipped with **safety buttons**.



### NUMBERED AND PADDED STRAPS

The straps are **numbered** to facilitate the correct closing sequence. Internally they have **MTP padding** and are **seamless**, for more comfort on the skin.

### MEDIAL CONDYLAR PAD

Anatomically designed and pre-assembled, to provide greater stability.

### TIBIAL PADDING

The tibia is positioned between two additional lateral shims, to enhance **the comfort and the stability** of the brace.

### SET OF CONDYLAR PADS

A practical set of pads, of **different thicknesses**, is included, suitable for any need.

### CLICKCLACK JOINT

It adjusts the flexion-extension by means of wedges, **without the need of any tool** (screwdriver or allen wrench). Polycentric and with a very low profile, you open it, adjust it, and close it.  
Extension: 0°, 10°, 20°, 30°, 40°,  
Flexion 45°, 60°, 75°, 90°  
Blocks at 0°, 10°, 20°, 30°.



### THERMOFORMED PADS

**Internally** (where in contact with the skin) they are covered with a hypoallergenic and non-slip fabric, which has a modern texture design **externally**. **Fastens** with low-profile snaps that are easily removed to clean the brace.



Never before  
has a knee brace  
been so easy to wear