



mediroyal 


Orthoses are our passion

Mediroyal Hook Clips and Orthotic Straps



Mediroyal Hook Clips and Orthotic Straps

Mediroyal has developed a product line with modular hook clips and strap material that can be easily integrated and adapted into orthotics.

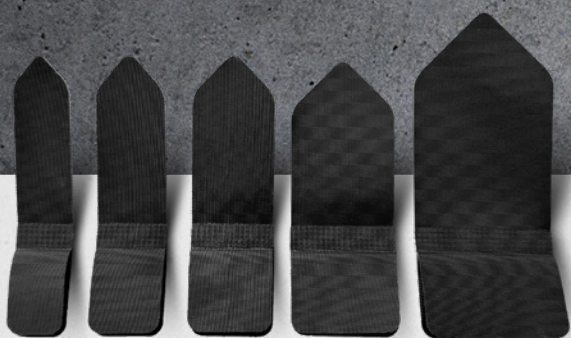

Mediroyal Hook Clips

The hook clips are made from an injection molded low profile hook with a strong grip that works both on conventional loop material as well as low profile OP textiles. The hook surface overlap is 30 mm in all widths.



Art.no	Width	Color	Quantity
EU445201	20 mm	Black	10 pcs/package
EU445202	25 mm	Black	10 pcs/package
EU445203	30 mm	Black	10 pcs/package
EU445204	38 mm	Black	10 pcs/package
EU445205	50 mm	Black	10 pcs/package

Mediroyal C-Loop Strap

C-Loop is a non-elastic strap based on a conventional loop strap that is lined on the backside and edges with a loop textile. The C-Loop strap is made like a tunnel, allowing you to insert reinforcements into the strap when riveting or sewing the strap to make the attachment point more stable.

Art.no	Width	Color	Roll
EU445212	25 mm	Black	15 m
EU445213	30 mm	Black	15 m
EU445214	38 mm	Black	15 m
EU445215	50 mm	Black	15 m





Mediroyal X-Strap

The X-Strap is developed by Mediroyal and features double sided OP loop material that has been welded together providing a stable strap with 2% reduced elasticity. This makes the strap useful in contoured orthotics where the straps need to run over complex shapes. The strap can be cut without fraying. The X-Strap is available in 10 m rolls in two widths.

Art.no	Width	Color	Roll
EU445240	40 mm	Black	10 m
EU445250	50 mm	Black	10 m

Mediroyal Double Soft Elastic

Double Soft Elastic is a thin, semi-elastic strap material that can be used secure and close smaller orthotic devices like foot parts on DAFO:s or resting hand splints. Double Soft Elastic conforms to the anatomy and can be attached onto the self-adhesive low-profile hook like EU445060 and EU445160. The Double Soft Elastic is available in 3,1 m rolls in two widths.



Art.no	Width	Color	Roll
EU450080	80 mm	Black	3,1 m
EU450110	110 mm	Black	3,1 m




Mediroyal Hook

The NRX® low profile hook has been specially developed to be used with our materials. They are injection molded with excellent grip strength and high durability. Available as self-adhesive, double laminated with hook on one side and loop on the top side and double-sided hook with standard hook on one side and strong hook on the other side.

The hook products are available in 25 m rolls

Art.no	Description	Width	Color
EU445060	Self-adhesive hook, standard	50 mm	Transparent
EU445160	Self-adhesive hook, strong	50 mm	Black
EU445061	Hook with loop top, standard	50 mm	Black
EU445161	Hook with loop top, strong	50 mm	Black
EU445260	Double sided hook, standard/strong	50 mm	Transparent/Black

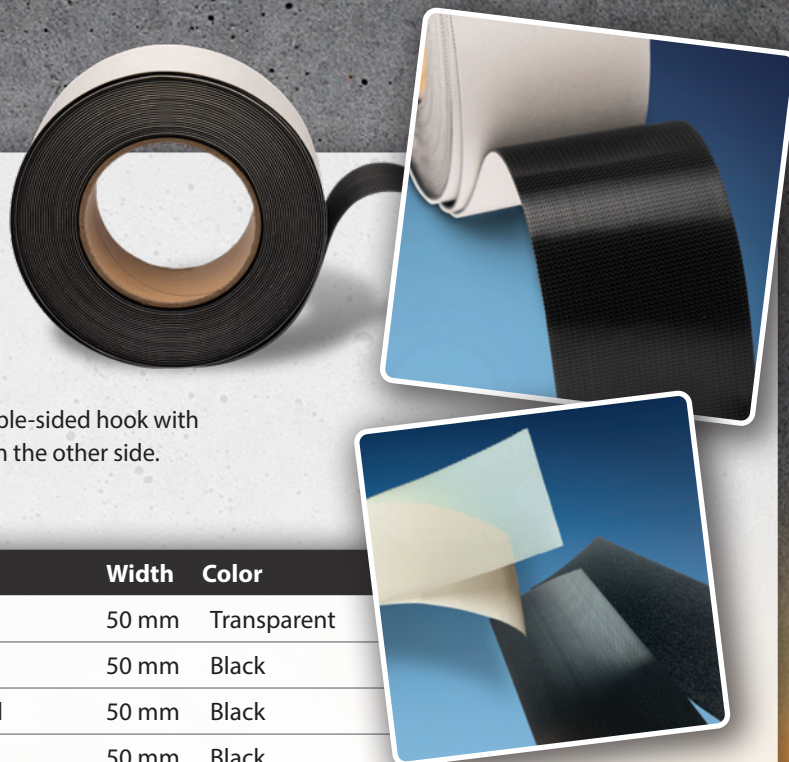

Orthoses are our passion

mediroyal 

Mediroyal Nordic AB
Box 7052 · Staffans väg 6B · SE-192 07 Sollentuna, Sweden
Tel 08-506 766 00 · info@mediroyal.se · www.mediroyal.se

Mediroyal Asia Pacific Co. Ltd.
Unit A10, 20F., No. 282, Shizheng N. 2nd Rd., Xitun Dist., Taichung City 407, Taiwan (R.O.C.)
+886-4-2254-8069 · sales.mediroyal@gmail.com · www.mediroyal.se/en